

SAFETY YOU CAN TRUST

ADVANCED TESTING

Helix test done as per IEC 13060 standard

CERTIFIED SAFETY

- Pressure and Temperature distribution test in side chamber through E-bro Logger
- Safety tests are done as per IEC 60601-1 Standard

TOUCH SCREEN CONTROL

Easy operation with accurate cycle monitoring



No. 40, 2nd Main Road, Pete Channappa Indl. Estate, Kamakshipalya,
Magadi Main Road, Bangalore 560 079.

 **80500 82373**

 Off.: 80500 82373 / 90710 87312

 sales@confidentelectronics.com

PHOS Class B Autoclave

Smart Touch Screen Operation
Modern Ergonomic Design
Robust Metal Cabinet



FEATURES ▼

Chamber & Materials


-  SS 316 Medical Grade Steel Chamber
-  Aluminum Instrument Tray
-  European Class B Standard with 4 Times Vacuum Drying

Safety Standards

-  Over Heat Control
-  Over Pressure Control
-  Safety Door Locking System
-  Calibrated Safety Valve

Advanced System

-  ARM Based Advanced Microcontroller System
-  User-Friendly Touch Screen Control
-  Inbuilt Steam Generator
-  Thermal Printer for Data Record

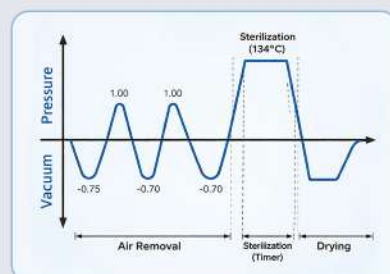

Advanced Sterilization Smart Control. Maximum Safety.



INBUILT THERMAL PRINTER
Automatic cycle data printout

PRINT DETAILS
Cycle Information

- Time • Temperature
- Pressure • Status



BASPA 23

TOUCH SCREEN DISPLAY & CONTROL PANEL

Display shows:

- Sterilization Mode
- Time
- Dry Time
- Vacuum Time
- Temperature
- Pressure
- Cycle Count



SECURE DOOR & HIGH-PERFORMANCE CHAMBER

Ergonomic front-loading door with secure safety locking and interlock system, heat and pressure-resistant sealing, and a high-grade SS 316 stainless steel chamber featuring uniform steam distribution, removable aluminium trays, and a smooth, corrosion-resistant interior for efficient sterilization and easy maintenance.

PRE PROGRAMMED STERILIZATION TABLE

PROGRAMME MODE	STERILIZATION TIME	DRY TIME	VACUUMING TIME	MAXIMUM LOAD	TOTAL TIME	CYCLE TYPE
BOWIE & DICK TEST	4 Minutes	4 Minutes	20 Minutes	Test pack	28 Minutes	Test Cycle
VACUUM TEST	-----	-----	15 Minutes	Empty	15 Minutes	Test Cycle
121 FAST	18 Minutes	5 Minutes	15 Minutes	5Kg	28 Minutes	S Cycle
121 WRAPPED SOLIDS	18 Minutes	10 Minutes	15 Minutes	2.5Kg	43 Minutes	S Cycle
121 POROUS	18 Minutes	15 Minutes	15 Minutes	2.5Kg	48 Minutes	B Cycle
121 HOLLOW WRAPPED	18 Minutes	10 Minutes	15 Minutes	2.5Kg	45 Minutes	B Cycle
134 FAST	5 Minutes	5 Minutes	15 Minutes	4Kg	25 Minutes	S Cycle
134 WRAPPED SOLIDS	5 Minutes	10 Minutes	15 Minutes	2.5Kg	30 Minutes	S Cycle
134 POROUS	5 Minutes	15 Minutes	15 Minutes	2.5Kg	35 Minutes	B Cycle
134 HOLLOW OPEN	5 Minutes	5 Minutes	15 Minutes	4Kg	26 Minutes	B Cycle
PRION (134°C)	20 Minutes	15 Minutes	15 Minutes	2.5Kg	50 Minutes	B Cycle

TECHNICAL SPECIFICATION (B CLASS AUTOCLAVES – 16L & 23L)

SPECIFICATION	CHAMBER CAPACITIES	
	16LTR	23LTR
OVERALL PRODUCT WIDTH	500 mm	500 mm
OVERALL PRODUCT HEIGHT	440 mm	440 mm
OVERALL PRODUCT LENGTH	610 mm	780 mm
CHAMBER DIAMETER	236 mm	236 mm
CHAMBER LENGTH	370 mm	540 mm
BASE BUSH (W × W)	430×490 mm	580×490 mm
MAX LOAD (NON-VACUUM)	4 Kg	5 Kg
INBUILT WATER TANK CAPACITY	2.5 Ltr	2.5 Ltr
NO. CYCLES PER FULL TANK	6-8	4-5
INSTRUMENT RACK PLATE SIZE	11.5×6.7×0.6 inches	18.3×6.7×0.6 inches
CHAMBER OPERATING SPECIFICATION		
OPERATING TEMPERATURE	0–148°C	0–148°C
OPERATING PRESSURE	0–2.45 bars	0–2.45 bars
ELECTRICAL SPECIFICATIONS		
VOLTAGE / WATTAGE	230V / 2800W (16L)	230V / 3050W (23L)
BAND HEATER ISOLATION	3.5 KV	3.5 KV
FREQUENCY	50–60 Hz	50–60 Hz
NOISE LEVEL (1M DISTANCE)	53.6 dBA	53.6 dBA
NOISE LEVEL (FRONT)	62.2 dBA	62.2 dBA
CHAMBER COMPONENT MATERIAL		
VESSEL	Stainless Steel 316 – S15	Stainless Steel 316 – S15
LID	Stainless Steel 316 – S15	Stainless Steel 316 – S15
TRAY	Aluminium 6061	Aluminium 6061